## **MEDICAL GRADE FABRIC MASKS**

Return with tracking BARCODE for Hospital: Hospital PROPERTY.


www.lovemendingsouls.org

1. Remove all stickers, place barcode in bag to return with masks. Fold fabric for double cutting. Cut straps out of same fabric, (around) 1/2" x 35". Each strap will be SINGLE ply so this makes 4 straps.

NOTE: If you leave the blue and white together while cutting, you will create four (4) straps from the two (2) cuts made from the wrap.

- 2. Fold fabric for double cutting. Keep white and blue wrap together and each pattern cut will make two (2) of the outline below, keep together, clips ok, NO PINNING! Each cut produces one full mask of left and right pieces. Best to cut out 3-5 of the pattern, and to place side-by-side, but flip them 180 degrees so it almost makes a straight line. This minimizes the waste.
- 3. Layer two (2) cut pieces together with <u>BLUE side out, leaving raw edge on outside.</u>
- 4. Sew with a straight stitch along the curved edge (1/4" seam allowance) with Polyester thread.

NOTE: BE SURE TO SECURELY BACKSTITCH THE BEGINNING AND ENDS OF ALL STITCH LINES

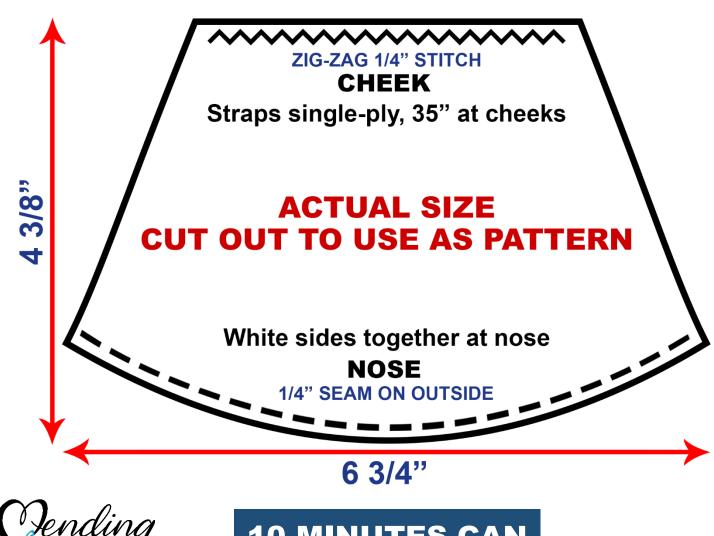

Sending Souls


10 MINUTES CAN SAVE A LIFE!

(520) 437-6007

## **COMFORT-FIT LAYOUT OF AUTOCLAVABLE MASKS**

From 24-inch square, 2-ply Blue-Blue Sheets Yields 8 masks with 4 ties each and 1 mask without ties Connect two 18" straps to make a 36" strap Attach this new 36" strap to each side of mask




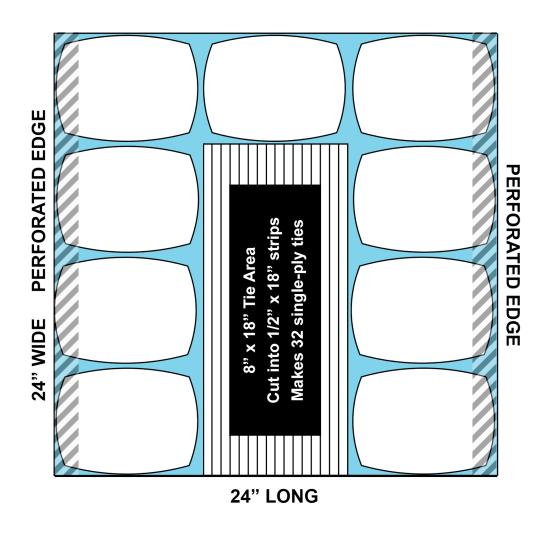

## **COMFORT-FIT LAYOUT OF AUTOCLAVABLE MASKS**

From 36-inch square, 2-ply Sheets (folded in half) Yields 18 masks and 24 pairs of tie straps 2 single-ply ties = 1 pair



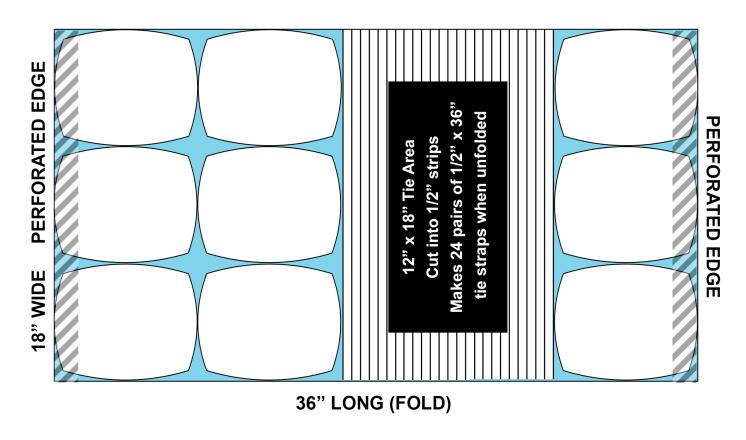

**NOTE:** If sheet is not full 36"x36", then tie strap area can be as small as 9" x 18" at the fold for 18 sets of 1/2" x 36" straps

## **COMFORT-FIT LAYOUT OF AUTOCLAVABLE MASKS**

From 48-inch square, 2-ply Blue-White Sheets
Fold at middle to cut
Yields 32 masks with ties and 4 masks without tie straps



